Unlimited
Combinations

Betriebsanleitung
Operating instructions
Notice d'instructions

Hei-SHAKE Reax 20

 **heidolph**
research made easy

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Seite 4 – 18

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Datum _____

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Firmenstempel _____

Translation of the original instructions

Page 22 – 36

Traduction de la notice originale

Page 40 – 54

Zertifikate/Certifications

Contents

Introduction

About this document..... 22

 Typographic conventions..... 22

 Copyright protection..... 22

 Safety instructions in the official languages of the European Union 22

General notes

Directives applied, product certification..... 23

Copyright notice 23

California Residents 23

Residual risk 23

Intended use 23

Compliant use 23

Reasonably foreseeable misuse 24

Transportation 24

Storage 24

Acclimatization 24

Permissible ambient conditions 24

Safety

General safety instructions..... 24

Electrical safety..... 25

Data security 25

Operational safety 25

Occupational safety 25

Personal protective equipment (PPE)..... 26

Environmental protection 26

Biohazard..... 26

Other regulations..... 26

Device description

Mechanical design..... 27

Commissioning

Set up the device.....	27
Power supply	27

Operation

Switching the device on/off	28
Equipment	28
Rotational mode.....	29

Troubleshooting

Troubleshooting.....	30
----------------------	----

Appendix

Technical specifications	31
Accessories	32
Attachments for bottles of 500 ml and 1,000 ml	32
Device service	33
General cleaning instructions.....	33
Replace the device fuse.....	33
Repairs - Return of equipment.....	34
Maintenance.....	34
Disposal.....	34
Warranty statement.....	34
Contact details	35
Certificate of decontamination.....	36

About this document

This operating manual describes the features and operation of vortexers, type Hei-SHAKE Reax 20. The operating instructions are an integral part of the described device!

Typographic conventions

This document uses the following symbols, signal words, and highlights:

	Description
	Warning symbols in combination with a signal word indicate dangers: DANGER Indicates an immediate dangerous situation. Failure to respect the indications will result in death or serious injury. WARNING Indicates a potential danger. Failure to respect the indications will result in serious injuries. CAUTION Indicates a potential hazard which, if not avoided, damage to property and minor to moderate injuries can occur.
	Mandatory signs are used to indicate important and useful information on handling a product. This information is used to ensure operational safety and to maintain the value of the product.
[...]	Square brackets identify labels of control elements on the device as well as labels and entries in software masks and user interfaces.
→	The arrow indicates specific instructions to be followed to ensure operational safety when handling the product.

Copyright protection

This document is protected by copyright and is intended for use by the purchaser of the product only.

No transfer to third parties, reproduction in any form, including excerpts, and by any means, as well as utilization and/or disclosure of the contents is permitted without the prior written consent of Heidolph Scientific Products GmbH. Any violation is subject to compensation for damage.

Safety instructions in the official languages of the European Union

A summary of all safety instructions in this document, in all official languages of the European Union, can be found in our Safety Guide for the Hei-PLATE product group (Ref. 01-005-006-68). This document is available for download on our homepage in the most up-to-date version.

Directives applied, product certification

CE marking

The device meets all requirements of the following directives:

- European Machinery Directive, 2006/42/EC
- EMC Directive, 2004/108/EC

Copyright notice

The software implemented in this product is protected by copyright laws. The rights holder is Heidolph Scientific Products GmbH, Walpersdorfer Straße 12, 91126 Schwabach/Germany. Any open source software components in this software are excluded from our copyright. Further information is available in the service area on our website www.heidolph.com.

Violations of copyright (e.g. unauthorized use or modification of the software) may result in civil claims (e.g. omission, compensation for damages) and/or entail criminal penalties within the meaning of the legislation of the destination country.

California Residents

Important information for California residents regarding Prop 65. Please visit www.P65Warnings.ca.gov for more information.

Residual risk

The device has been designed and manufactured in accordance with the state-of-the-art standards at the time of development and the recognized safety regulations. During mounting and use, as well as during maintenance, repair and cleaning work, there are nevertheless certain residual risks associated with the described device.

These are identified and described at the appropriate points in this document.

Intended use

The described platform shaker in this document has been specially developed for the following tasks:

- Shaking
- Mixing
- Emulsifying
- Suspending
- Separating

The areas of application of the described device include chemical, biological and environmental-analytical laboratory and research applications, applications of basic research and comparable facilities.

Due to its design, the device in its delivery condition may only be used in analytical processes or in laboratory-like conditions in the food, cosmetics, and pharmaceutical industries as well as other comparable industries that manufacture products intended for consumption by humans or animals, or for use on humans or animals.

Any other use of this device is not considered as intended!

Compliant use

The user is generally responsible for evaluating the conformity of his application and, if necessary, for taking additional measures.

Reasonably foreseeable misuse

Additional measures may be necessary for use under conditions or for purposes deviating from the intended use, and/or specific guidelines and safety regulations must be observed. Corresponding requirements must be evaluated and implemented by the operator in each individual case.

Compliance with and implementation of all relevant directives and safety measures for the respective field of application is the sole responsibility of the operator.

All risks resulting from improper use are borne solely by the operator.

The device may only be operated by authorized and instructed personnel. Training and qualification of the operating personnel as well as ensuring that the device is handled responsibly are the sole responsibility of the operator!

Transportation

During transportation, avoid severe shocks and mechanical stresses that can cause damage to the device.

Keep the original packaging in a dry and protected place for later use.

Storage

Always store the device in its original packaging. To protect against damage and excessive material aging, store the device in an environment that is as dry, temperature-stable and dust-free as possible.

Acclimatization

After each transportation and after storage under critical climatic conditions (e.g. high temperature difference between inside and outside), allow the device to acclimatize at room temperature for at least two hours to prevent possible damage due to condensation before commissioning it in the place of use. If necessary, extend the acclimatization phase if the temperature differences are very high.

Make all supply connections (power supply, tubing) only after the device has been acclimatized!

Permissible ambient conditions

The device is designed for indoor use only. The device is **NOT** suitable for outdoor use!

The device is **NOT** suitable for use in potentially explosive areas!

When used in corrosive atmospheres, the service life of the device may be reduced depending on the concentration, duration and frequency of exposure.

General safety instructions

Before commissioning and using the device, familiarize yourself with all the safety regulations and occupational safety guidelines applicable at the place of use and observe them at all times.

Only operate the device if it is in faultless technical condition. In particular, make sure that the device is free of visible damage.

If there is missing or misleading information on the device or regarding occupational safety, contact the responsible safety specialist or our technical service.

Only use the device in accordance with the regulations on intended use („Intended use“ on page 23).

Electrical safety

Before connecting the device to the power supply, ensure that the voltage indicated on the rating plate matches the specifications of the local power utility company.

Ensure that the power supply circuit provided is protected by means of a residual-current device (RCD).

Do not make any changes to the hard-wired three-pole power supply (PH, N, GND).

Always replace defective device fuses with fuses of the same type with identical performance features.

Prior to use, check that the device and the power supply cord are free of visible damage.

Have repairs and/or maintenance work on the device carried out exclusively by an authorized electrician or by the technical service department of Heidolph Scientific Products GmbH.

Always switch off and disconnect the device from the power supply, preventing reconnection, before carrying out maintenance work, cleaning, or repairs.

Data security

The user is responsible for ensuring data security when transferring data between the described device and other devices.

Only use secure networks for the data transfer and avoid use of critical infrastructure.

Only use high-quality shielded data cables for the data transfer.

For data transfer via a USB B connection, an industrial standard USB hub should be preferably used to ensure the most stable connection possible.

Operational safety

Operate the device under a closed ventilated fume hood when working with potentially hazardous substances (see EN 14175 and DIN 12924).

Do not make any unauthorized changes or modifications to the device!

Only use genuine spare parts and accessories, or those expressly approved by the manufacturer!

Rectify malfunctions or faults on the device immediately.

Switch off and disconnect the device from the power supply, preventing reconnection, if it is not possible to eliminate the malfunction or rectify the fault immediately.

Observe all other applicable regulations such as laboratory and workplace guidelines, recognized safety technology rules and special local regulations.

Occupational safety

Always use the prescribed personal protective equipment (PPE) such as protective clothing, safety goggles, protective gloves, safety shoes, etc.

Do not operate any other devices in the immediate vicinity of the device ...

- which can generate electromagnetic fields in the frequency range between 9×10^3 Hz to 3×10^{11} Hz,
- which generate emission or radiation sources in the frequency range 3×10^{11} Hz to 3×10^{15} Hz (in the optical spectral range wavelengths from 1,000 μm to 0,1 μm),
- which generate ultrasonic or ionizing waves.

Do not process any substances that could release energy in an uncontrolled manner (e.g. self-ignition).

Do not process substances in which the energy input through mixing poses a danger.

Wipe off any liquid that may have spilled on the device immediately.

Personal protective equipment (PPE)

The operating company must determine and provide the necessary PPE, depending on the respective application and the media and chemicals used.

The corresponding instruction of the personnel is solely within the operating company's responsibility.

Environmental protection

When processing environmentally hazardous substances, take appropriate measures to avoid hazards to the environment.

The evaluation of corresponding measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Biohazard

When processing biohazardous substances, take appropriate measures to prevent hazards to persons and the environment, including:

- Instruction of the personnel regarding the necessary safety measures.
- Provision of personal protective equipment (PPE) and instruction of the personnel in its use.
- Marking of the device with a biohazard warning symbol.

The evaluation of corresponding measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Other regulations

In addition to the notes and instructions in this document, observe all other applicable regulations such as laboratory and workplace guidelines, hazardous substances ordinances, recognized rules of safety engineering and occupational medicine as well as particular local regulations!

Noncompliance will invalidate any warranty claims against Heidolph Scientific Products GmbH.

The operator is generally liable for all damage resulting from unauthorized changes or modifications to the device, from the use of unapproved or non-genuine spare parts and accessories, or from disregarding the safety instructions and hazard warnings or the manufacturer's instructions!

Mechanical design

Front view

The figure shows an example of the device variant Hei-SHAKE Reax 20/8:

- 1 Clamping lever (upper side and underside)
- 2 Emergency stop button
- 3 On/Off switch
- 4 Speed controller

Set up the device

CAUTION

Improper installation of the device may result in direct and indirect damage to property due to falling and/or spilling of fluids! Observe the instructions for the correct positioning of the device. Observe the generally valid rules for the safe installation of work equipment.

The surface must be even, non-slip, proper and sufficiently stable.

The device must be freely accessible at all times.

The professional mounting and positioning of the device including the supplied accessories is within the sole responsibility of the operator! Heidolph Scientific Products GmbH accepts no responsibility for direct and/or indirect personal injury or damage to property resulting from failure to observe the instructions for the correct mounting and positioning of the device.

Power supply

DANGER

Please note the safety instructions in section „Electrical safety“ on page 25

The hard-wired power supply leaves from the rear of the device. Make sure that the device is switched off (main switch in position off).

Switching the device on/off

Use the main switch [Power] at the front to turn the device on or off.

Equipment

In delivery status, the device is equipped with pockets for bottles of 2,000 ml. You will need special pockets for bottles of 500 ml and 1,000 ml. These are available as optional accessories (for options and mounting, see „Attachments for bottles of 500 ml and 1,000 ml“ on page 32).

To insert the bottles, you first need to move the rotor in a suitable position:

- Adjust the speed controller to the lowest rotation speed [1].
- Switch on the device.
- Switch off the device again as soon as the rotor is in a suitable position.

- Loosen one clamping lever [1] (turn counterclockwise around 90°).


- Pull up the knurled handle of the retaining bolt on the loosened clamping lever [2] and insert the bottle.
 - Make sure to place the pressure plate centrally in relation to the bottleneck.

- Press down the retaining bolt to fix the bottle with the pressure plate:
- Lock the loosened clamping lever again (turn clockwise around 90°).

- If necessary, the clamping lever position can be adjusted: to do so, pull the clamping lever outwards, turn the clamping lever in a suitable position and loosen the clamping lever again.
 - Repeat the process, if necessary, until being able to clamp the bottle correctly.

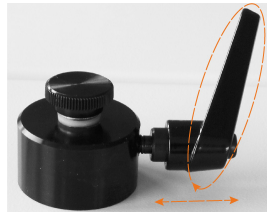Rotational mode

WARNING

Always place individual vessels in the center. Always ensure that the weight load is evenly distributed.

Use suitable pockets for the used vessels. Further information on available accessories can be found on our website at www.heidolph.com!

Before starting the rotational movement, make sure that all bottles are correctly closed and securely clamped.

Reduce the rotational speed if the device starts to vibrate during operation.

Always use the necessary and appropriate personal protective equipment!

Never touch rotating parts with your hands. Do not bend over rotating parts.

Wear tight-fitting work clothes, remove jewelry before starting work, tie long hair under a hair net or under a suitable head covering

- Adjust the speed controller to the desired rotation speed (adjustment range: 1 – 16 rpm).
- Switch on the device.
 - The rotational movement starts with the lowest rotation speed and will gradually be increased to the set speed.
 - Observe the stability of the device during its start-up and until reaching the set rotation speed.
 - The rotation speed can be readjusted continuously at any time.

In case of a power interruption, the device has to be switched on again with the on/off switch when the mains return.

Troubleshooting

Failure	Cause	Possible remedy
Rotational movement does not start	No mains voltage	Check if the on/off switch is in position [on]
	Power supply cord defective	Check the power supply cord for damages/correct connection
	On/off switch in position [off]	Check building fuses Check device fuses
Rotational movement stops unexpectedly	Motor overheat protection has tripped	Wait about 20 minutes and reduce the load
	Electrical fault (no engine noise)	Contact Heidolph Sales or Service
	Mechanical defect (engine noise audible)	

In case of recurring errors please contact the responsible sales department or our technical service (see section „Contact details“ on page 35).

Technical specifications

Hei-SHAKE Reax 20	Reax 20/4	Reax 20/8	Reax 20/12
Dimensions (W × H × D) in mm	490 × 520 × 465	770 × 520 × 465	1050 × 520 × 465
Weight (kg)	23	28	33
Permissible load	max. 30 kg		
Bottle dimensions	Maximum bottle diameter: 136 mm Maximum bottle height: 270 mm		
Motion	Overhead, rotating		
Speed range/speed display	1 – 16 rpm/scale		
Speed control	electronic		
Drive	Asynchronous motor		
Motor protection	Overheat protection, self-resetting		
Protection class	IP21 (acc. to EN 60529)		
Protection class	I ⚡ (acc. to EN 61140)		
Acoustic pressure	< 45 dB (A), in accordance with IEC 61010		
Rated voltage	230 V, 50/60 Hz 115 V, 60 Hz		
Fuse	T3, 15 A (230 V) T6, 25 A (115 V)		
Power input	280 W		
Operating and storage temperature	5 °C – 31 °C at up to 80 % rel. humidity 32 °C – 40 °C at up to 50 % rel. humidity (decreasing linearly)		
Installation altitude	up to 2,000 m asl		

Scope of delivery

Component	Quantity	Product no.
Hei-SHAKE Reax 20/4	1	541-20004-00
Hei-SHAKE Reax 20/8	1	541-20008-00
Hei-SHAKE Reax 20/12	1	541-20012-00
Power supply cord	1	country specific
Warranty registration	1	01-006-002-78

Accessories

Component	Quantity	Product no.
Device fuse 230 V	1	14-002-015-46
Device fuse 115 V	1	14-002-015-49
Attachment for bottles of 500 ml Bottle height 162 – 202 mm	1	549-27000-00
Attachment for bottles of 1,000 ml Bottle height 208 – 248 mm	1	549-26000-00
Tension plate for sealing caps Ø 94 mm	1	11-001-001-81

Attachments for bottles of 500 ml and 1,000 ml

Example: attachment for bottles of 1,000 ml, consisting of 4 × spacer bolt and 1 × hex wrench.

- Tighten the clamping lever.
- Loosen the pressure plate using the hex wrench.

- Screw the spacer bolt with the retaining bolt (note correct position of the tension spring).

- Screw the pressure plate with the extended retaining bolt.

Further information on available accessories can be found on our website at www.heidolph.com!

Device service

When carrying out service work on the device (cleaning, maintenance, repair), observe the general instructions and safety information described in this section.

DANGER

Electric shock

Live components are installed inside the device. When opening the device, there is a risk of touching live components.

Switch the device's main switch off and disconnect it from the power supply before carrying out maintenance work, cleaning, or repairs.

Penetrating liquid poses the risk of an electric shock.

When cleaning, avoid the penetration of liquids.

General cleaning instructions

Wipe all surfaces and the control panel with a damp cloth if necessary. Persistent contamination can be removed with mild soapy water.

CAUTION

Property damages

Improper cleaning can damage the surfaces of the device. Penetrating liquid can damage the electronic components inside the device.

Clean the device's surfaces with a soft, lint-free and only slightly moistened cloth.

Never use any aggressive or abrasive cleaning agents or aids.

Replace the device fuse

The fuse holder for the two device fuses is located on the rear of the device below the IEC appliance inlet.

DANGER

Before replacing the fuses, switch off the device and disconnect it from the power supply.

Always replace the two device fuses in pairs with original manufacturer fuses. Further information on available accessories can be found on our website at www.heidelph.com!

After fuse replacement, check the device for a safe condition according to IEC 61010-1.

Repairs - Return of equipment

Repairs to the device may only be carried out by authorized experts! Unauthorized repairs during the warranty period will result in the loss of the warranty claim. Regardless of the warranty claim, the owner is solely liable for damage caused by unauthorized repairs.

- In case of repair and before returning your device, contact our technical service at the following e-mail address:
 - service@heidolph.de.
- In your message, please provide us with the following information in addition to an error description:
 - Item number
 - Serial no.

The required data can be found on the rating plate of the device.

A service representative will contact you as soon as possible to agree on the next steps.

Include the completed certificate of decontamination with every device return, see „Certificate of decontamination“ on page 44.

Maintenance

The device contains no user-serviceable components. If necessary (in the event of abnormal operating behavior such as excessive noise or heat generation, for example) contact our technical service.

Disposal

- When disposing of the device, observe the provisions of the WEEE Directive 2012/19/EU and its transposition into national law in the country of use.
- When disposing of portable batteries, observe the provisions of the European Battery Directive 2013/56/EU and their implementation in national law in the country of use.
- Check the device and all components for residues of substances that are hazardous to health, the environment and biohazardous before disposal.
- Properly remove and dispose of residues of substances that are hazardous to health, the environment and biohazardous!

Warranty statement

Heidolph Scientific Products GmbH provides a three-year warranty against material and manufacturing defects.

Glass and wear parts, transportation damage, and damage resulting from improper handling or non-intended use of the product are excluded from the warranty.

The warranty period for registered products begins on the date of purchase. Register the product with the enclosed warranty card or on our homepage www.heidolph.com.

For non-registered products, the warranty period begins with the date of the serial production (to be determined by the serial number).

In the event of material or manufacturing defects, the product will either be repaired or replaced free of charge within the warranty period.

Contact details

Heidolph Scientific Products GmbH

Technical service

Walpersdorfer Str. 12

D-91126 Schwabach/Germany

Email: service@heidolph.de

Representations

To find your local Heidolph distributor please visit www.heidolph.com

Certificate of decontamination

Enclose the Certificate of decontamination, duly completed, with your device return. Submissions without a Certificate of decontamination cannot be processed!

CERTIFICATE OF DECONTAMINATION IN CASE OF RETURNS

Please fill in the required fields.

Note: The sender must package the goods properly and appropriately for transport.

Heidolph Scientific Products GmbH

Walpersdorfer Straße 12
91126 Schwabach

Phone: +49 (0) 9122 9920-380

Fax: +49 (0) 9122 9920-19

E-Mail: service@heidolph.de

SENDER

Name _____	First name _____
Company/institution _____	Department _____
_____	Workgroup _____
Address _____	
ZC/City _____	
Country _____	Phone _____
Email _____	

DEVICE DETAILS

Article number _____	Serial no. _____
Ticket number _____	
Reason for sending in _____	

Has the device been cleaned, decontaminated/disinfected? **Yes** **No** (Please mark as applicable)

If yes, which measures were carried out?

Does this device pose a risk to people and/or the environment due to the processing of substances that are hazardous to health, the environment and/or are biohazardous? **Yes** **No** (Please mark as applicable)

If yes, with which substances did the device come into contact?

LEGALLY BINDING DECLARATION

The principal/consignor is aware that they are liable to the agent/consignee for losses or damage incurred due to incomplete and incorrect information.

_____	_____	_____
Date	Signature	Company stamp

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hereby declare, that the product (from serial number on 200066855) designated below is in compliance with the basic requirements of all applicable EU-directives stated below with regard to design, type of model sold and manufactured by us. This certificate will be invalid if the product is modified without the prior written consent and agreement of the manufacturer.

Reax 2	541-21001-XX		
Reax top	541-10000-XX		
Reax control	541-11000-XX		
Reax 20/4	541-20004-XX		
Reax 20/8	541-20008-XX		
Reax 20/12	541-20012-XX		
Titramax 100	544-11100-XX, 544-11200-XX		
Duomax 1030	543-32105-XX, 543-32205-XX, 543-32210-XX		
Polymax 1040	543-42105-XX, 543-42110-XX, 543-42205-XX, 543-42210-XX		
Polymax 2040	542-40005-XX, 542-40010-XX		
Multi Reax	545-10000-XX		
Vibramax 100	544-21200-XX	Titramax 101	544-11300-XX
Rotamax 120	544-41200-XX	Vibramax 110	544-31200-XX
Titramax 1000	544-12200-XX	Inkubator 1000	549-90010-XX
Unimax 1010	543-12310-XX	Promax 1020	543-22332-XX
Unimax 2010	542-10020-XX	Promax 2020	542-20020-XX
Hei SHAKE Orbital Core	546-11XXX-XX		

Maschinenrichtlinie / Machinery Directive 2006/42/EG
EMV-Richtlinie / Electromagnetic Compatibility Directive 2014/30/EU
Delegierte (RoHS-) Richtlinie / Delegated (RoHS) Directive 2015/863/EU
Angewandte (harmonisierte) Normen / (Harmonized) Standards applied:
EN ISO 12100:2010, EN 61326-1:2013, EN 61010-1:2010, EN IEC 63000:2018

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Jörg Ziel
Quality Manager

Declaration of Conformity

In accordance with UK Government guidance

This declaration is issued under the sole responsibility of the manufacturer,

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Hei SHAKE Orbital Core 546-11XXX-XX

Description:

Shaker & Mixer with different working directions

The object of the declaration described above is in conformity with the relevant UK Statutory Instruments (and their amendments):

2008 No. 1597

2008 2016 No. 1091

2012 No. 3032

The Supply of Machinery (Safety) Regulations

The Electromagnetic Compatibility Regulations 2016

The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012

and complies with the following technical standards :

EN ISO 12100:2010, EN 61326-1:2013, EN 61010-1:2010, EN IEC 63000:2018

UK Authorised Representative (for authorities only)

ProductIP (UK) Ltd.

8. Northumberland Av.

London WC2N 5BY

Signed for and on behalf of Heidolph Scientific Products GmbH

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With this declaration, we confirm (according to current knowledge) that all sold laboratory devices and components by Heidolph Scientific Products GmbH fulfill the requirements of the EU directive 2011/65/EU (RoHS) and the extension 2015/863. All devices are compatible with the requirements of the RoHS Directive for the following materials:

Max. 0,01% of weight in homogenous material for cadmium and max. 0,1% of the weight in homogeneous material for Lead, Mercury, Hexavalent chromium, Polybrominated biphenyls (PBB), Polybrominated diphenyl ethers (PBDE), Bis(2-ethylhexyl) phthalate (DEHP), Butyl benzyl phthalate (BBP), Dibutyl phthalate (DBP), Diisobutyl phthalate (DIBP).

In the case of individual assemblies, maximum concentrations maybe exceeded within the permissible exceptions of the Directive.

Schwabach, 04.September.2024

Andreas Hahn
Managing Director

Jörg Ziel
Quality Manager

China RoHS DECLARATION OF CONFORMITY

Heidolph Scientific Products GmbH has made reasonable efforts to ensure that hazardous materials and substances may not be used in its products.

In order to determine the concentration of hazardous substances in all homogeneous materials of the subassemblies, a "Product Conformity Assessment" (PCA) procedure was performed. As defined in GB/T 26572 the "Maximum Concentration Value" limits (MCV) apply to these restricted substances:

- Lead (Pb): 0.1%
- Mercury (Hg): 0.1%
- Cadmium (Cd): 0.01%
- Hexavalent chromium (Cr(VI)): 0.1%
- Polybrominated biphenyls (PBB): 0.1%
- Polybrominated diphenyl ether (PBDE): 0.1%

Environmental Friendly Use Period (EFUP)

EFUP defines the period in years during which the hazardous substances contained in electrical and electronic products will not leak or mutate under normal operating conditions. During normal use by the user such electrical and electronic products will not result in serious environmental pollution, cause serious bodily injury or damage to the user's assets.

The Environmental Friendly Use Period for Heidolph Scientific Products GmbH products is 25 years.

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This table is created according to SJ/T 11364-2014

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Note: Table applies to all products. Some of the components or parts listed above may not be part of the enclosed product.

- O: 表示该有毒有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。
O: Indicates that the above mentioned hazardous substance contained in all homogeneous materials of the part is below the required limit as defined in GB/T 26572.
- X: 表示该有毒有害物质至少在该部件某一均质材料中的含量超出GB/T 26572规定的限量要求。
X: Indicates that the above mentioned hazardous substance contained in at least one of the homogeneous materials of this part is above the required limit as defined in GB/T 26572.

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Apart from the disclosures in the above table, the subassemblies are not intentionally manufactured or formulated with lead (Pb), mercury (Hg), cadmium (Cd), hexavalent chromium (CrVI), polybrominated biphenyls (PBB), and polybrominated diphenyl ethers (PBDE).

Products manufactured by Heidolph Scientific Products GmbH may enter into further devices or can be used together with other appliances.

With these products and appliances in particular, Heidolph Scientific Products GmbH will not take responsibility for the EFUP of those products and appliances.

Schwabach, 04.September.2024

Andreas Hahn
Managing Director

Jörg Ziel
Quality Manager

© Heidolph Scientific Products GmbH

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Subject to change without notice. The printed version of this document is not regularly updated. The latest issue of this document can be found by visiting our homepage.

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MATERIAL CONTENT DECLARATION FOR Heidolph Scientific Products GmbH PRODUCTS